



# Start to Finish Central Line Placement

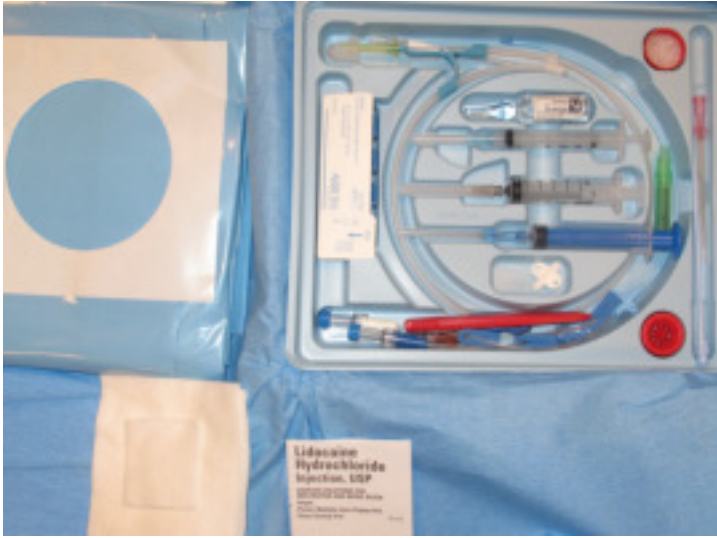


Triple Lumen Central Venous Catheter Kit

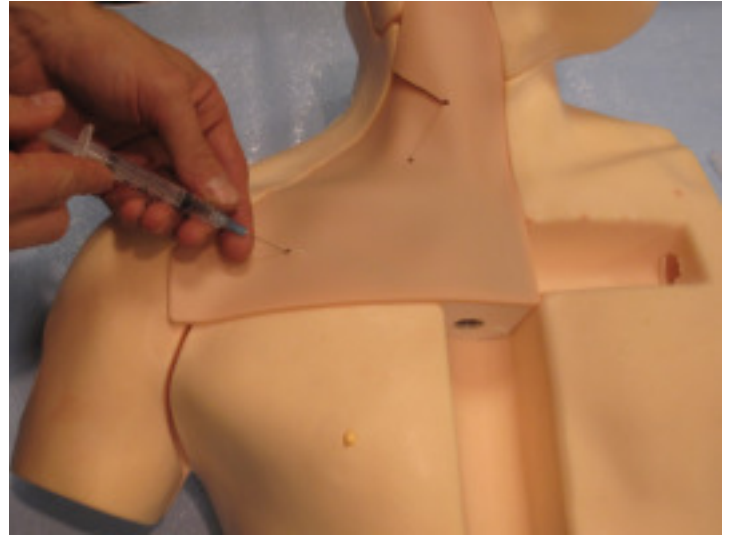


**Step 1:** Anesthetize the skin with lidocaine




**Step 2:** Advance introducer needle into vein. Flash of blood will NOT be seen

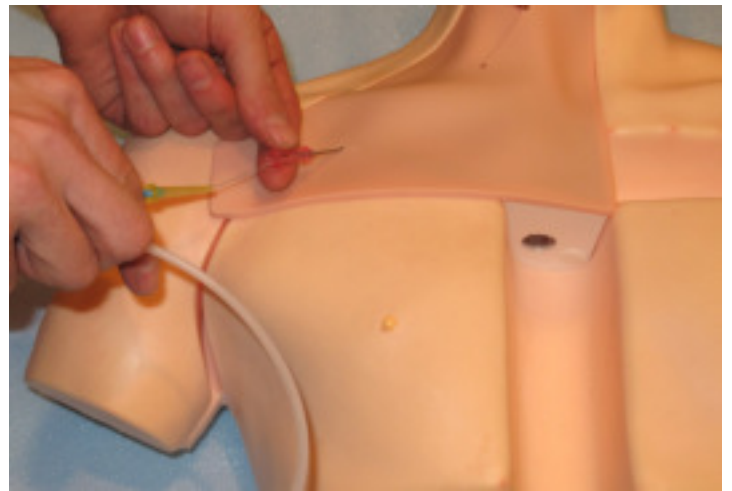

**Step 3:** Remove syringe and advance guidewire to 20 cm

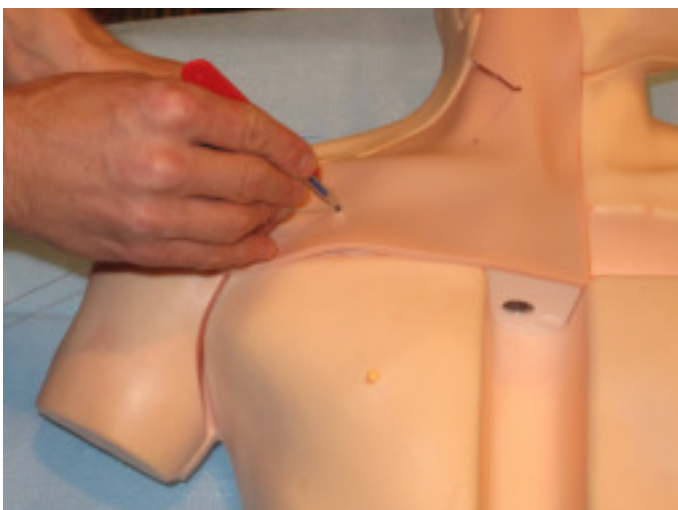


**Step 4:** Nick the skin with a scalpel over the wire

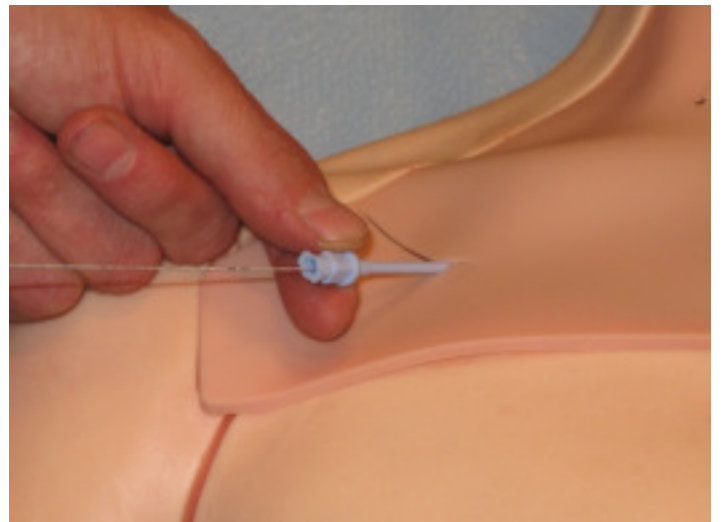


**Step 5:** Advance dilator over guidewire with slow twisting motion


# Start to Finish Central Line Placement

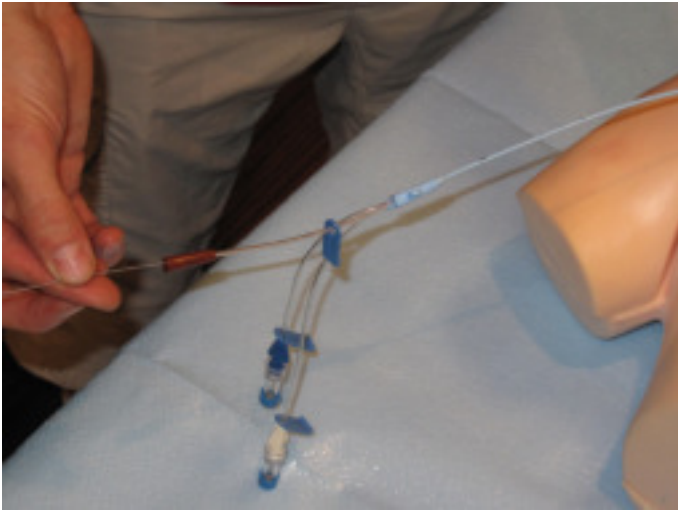

**Step 6:** Remove dilator and advance catheter over guidewire maintaining control of the guidewire at all times

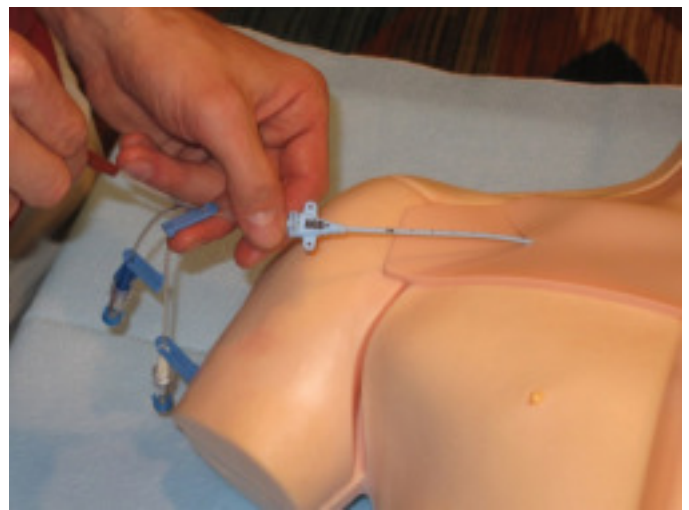

**Step 7:** Remove guidewire from brown port of catheter

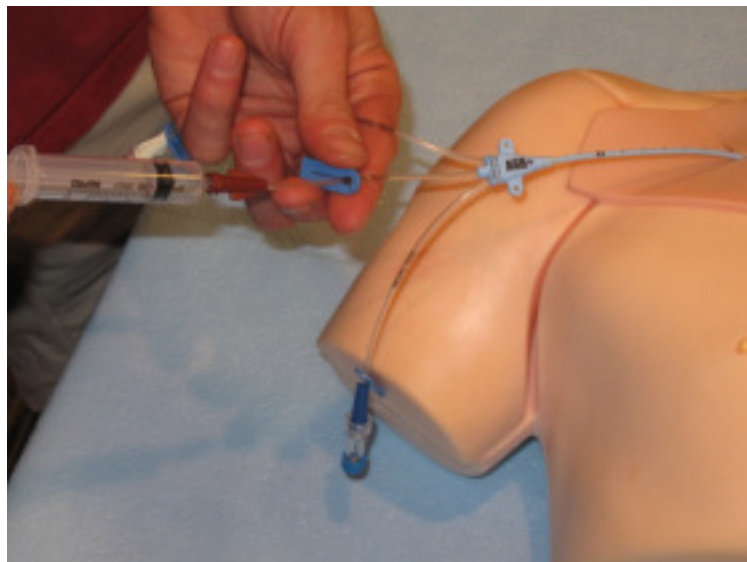

**Step 8:** Attach syringe, unclamp, aspirate then flush each catheter port